500 XLE comfort

The reclining treatment chair designed for a particularly high level of comfort and two-tone design.

Order No. 15000006 + 30370080

ADVANTAGES



- Optimized price-performance ratio
- Ample legroom and unobstructed access to the patient's head
- Space-saving design for flexible use, outstanding maneuverability even in the most confined of spaces
- Multi-functional manual keypad with 4 memory positions and auto-run function
- Automatic switch-off in the stand-by mode with reactivation button
- Sterile working conditions without interruption thanks to adjustment functions provided by the foot keypad
- Not location-bound thanks to battery operation
- The treatment chair featuring comfort upholstering is particularly well suited for longer surgical procedures (e.g. retina surgery)



- Easy maneuverability thanks to mechanic guidance operated via foot switch
- · 360-degree rotating, electrically conductive double casters for outstanding movement characteristics
- Stool & table combination no need for shifting or relocating the patients
- Service-friendly access for easy cleaning and maintenance
- User-friendly



- Easy mounting and dismounting, low access height
- Comfortable armrests for high patient safety, removable
- Patient-friendly comfort upholstery including integrated knee support cushion in a two-tone design























Memory manual switch with four freely storable memory positions and auto-run function.



Integrated knee support cushion, individually positionable to suit every patient's body height.



Two-tone comfort upholstery for even more comfortable sitting/reclining positions.



Zinc-primed chassis with anthracite structured finish. Scratch-resistant, hard-wearing, ABS covering in white.



Detachable push handles on the foot ends for easy maneuvering.



Telescopic stainless-steel accessories storage system.





Two-tone comfort upholstery (anthracite F6461197 / chrome F6461205 standard shade) for even more comfortable sitting/reclining positions.



Available in a wide variety of colors. Please refer to Page 53 for our complete range of colors.

500 XLE comfort

500 XLE comfort









ACCESSORIES

FLEXIBLE SURGICAL DRAPE SUPPORT ROD without oxygen connection
 No. 32050000
 with oxygen connection
 No. 32060000

2 HEAD PART BUMPER stainless steel No. 30280000

3 HAND CONTROL HOLDER, PVC No. 30390002

4 BACKUP BATTERY No. 30470000

5 SIDE PARTS, (right + left) swiveling

6 ADJUSTABLE WRIST WREST incl. mounting plate

MONITOR ARM stainless steel, to be swung,
 incl. I.V. Pole stainless steel Ø 25 mm and fixation
 plate monitor stand
 No. 30240014

PADDED ARM BOARD with universal joint,with wrist secure strapsNo. 33510011

IO I.V. POLE stainless steel ø 25/18 mm, one-hand height-adjustable 1070–1660 mm, 4 hooks, max. 2 kg 84.4 lbs/hook, glass and holder No. 30680000









No. 30140000






500 XLE comfort





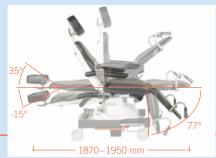



SPECIFICATIONS

Chassis

- Frame profile
 Outside length
 Outside width
 Width incl. armrest
 90 x 50 mm
 805 mm
 605 mm
 830 mm
- Zinc undercoat, elegant anthracite matte finish
- Brake operable from both sides and acting on two (2) wheels each; 2 wheels featuring independent directional wheel positioning
- Double rimmed casters ø 125 mm
- Detachable foot keypad (total up/down, head up/down)
- Emergency stop button
- Visible screws are stainless steel
- With equipotential bonding connection

Lifting Column

- Max. lifting force	300 mm
- Max. patient weight	250 kg
- Motors IPX4 with HAL sensoric, to be memorized	

Upper Frame

- 4-part support frame (head, back, seat, foot)
- Length 1,870–1,950 mm - Width 610 mm

- Access height
- Removeable upholstery
- Tundra skai or skai Pandoria Plus colors see color chart

May Load Angle

	Max. Load	Angle	Dimensions (i x w)
- Head part	20 kg	-15°/+35°	300 x 250 mm
- Back part	90 kg	0°/+85°	610 x 610 mm
- Seat part	90 kg	-5°/+35°	450 x 610 mm
- Leg part	50 kg	0°/+77°	510 x 610 mm
- Trendelenburg		-12°	
- Height adjustment	250 kg		from 580-880 mm
- Armrests	30 kg ea.	0°/+180°	(2 x) 470 x 75 mm

Electronic Control

- All components according to IPX4
- Manual keypad with 4 individually programmable memory positions, auto-run function, reset button
- Foot keypad

Battery unit

- Battery box with 2 batteries (lead-gel), pollution free, maintenance-free, cadmium-free 2 x 12 V 7.2 Ah
- As back-up-system exchangeable without tools

- Acoustic low voltage signal for discharge protection
- Universal battery charger Input: 100-240 V, Output: 24 V

Operating Time Per Battery Charge

- Measured by operating circle (Prepare-Operate-Recover) 30 0
- Quantity of operating circles depends on individual operation

Electromotive Drives

580 mm

- Head, back, seat, foot, total up/down
- Trendelenburg/Reverse Trendelenburg

Further Technical Details

- Chair has detachable, swiveling armrests made completely with PU foam
- All touchable metal parts are grounded
- Detachable seat for service friendly access to electric components

Max. Load	250 kg	
	•	
Weight	135 kg	

Subject to technical modifications.