

hana SSXT™

THE ULTIMATE SPECIALTY SURGERY EXTENSION FOR ANTERIOR APPROACH THA



MIZUHOSI®
BRINGING INNOVATION TO THE TABLE.™

hanaSSXT™

THE ULTIMATE SPECIALTY SURGERY EXT




JOIN THE SWITCH TO ANTERIOR APPROACH THA


The Anterior Approach to total hip arthroplasty (AA THA) has become the procedure of choice for many surgeons and patients alike. Since 2005, the hana® Orthopedic Table has been used to successfully position patients in thousands of these procedures. Knowing it as a vital instrument for the procedure, surgeons are reluctant to attempt AA THA without a hana® Table. The hanaSSXT™ now offers many of the features from the full hana® Table and is available to expand the capability of your general surgery table to successfully perform AA THA.

A PLATFORM DESIGNED FOR AA THA AND HIP ARTHROSCOPY

Easily mounted to most standard general surgery tables, the hanaSSXT™ provides precise patient positioning that enables the surgeon to replace the hip through a single anterior incision without detachment of muscle from the pelvis or the femur. It allows for the safe hyperextension, abduction, adduction and external rotation of the leg required for the AA THA procedure. This positioning capability is not available on any general surgery table today.

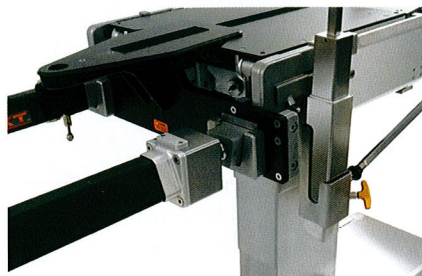


EXTENSION FOR ANTERIOR APPROACH THA

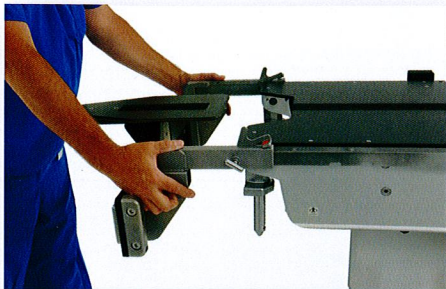


hanaSSXT™ FEATURES

- Fits most general surgery tables*
- Easy set-up and removal
- Exceptional maneuverability
- Extensive imaging capability
- Controlled femoral lift
- Precise patient positioning

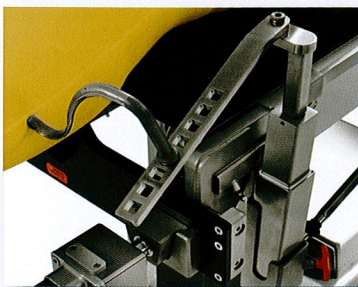

PROPRIETARY FEMORAL LIFT

- Allows precise femoral lift control
- Improves surgical access for femoral canal prep and femoral component placement

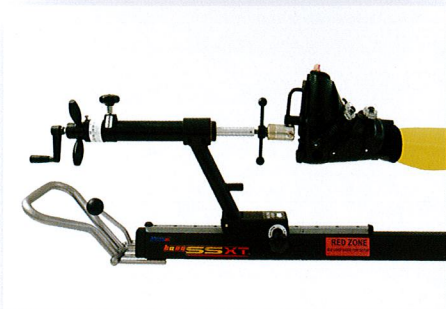


SPAR AND TRACTION BOOT ADVANTAGES

- Lightweight, carbon fiber spar construction
- Provides solid support of lower extremities
- Allows precise control and stability of lower extremity position, manipulation and traction
- Spars include Spherical Spatial Positioning System SSPS™



EXTENSIVE IMAGING CAPABILITY

- Radiolucent leg spars provide uninterrupted imaging
- Surgeon can intraoperatively confirm implant placement and accuracy of leg length
- Unrestricted C-Arm access to lower extremities

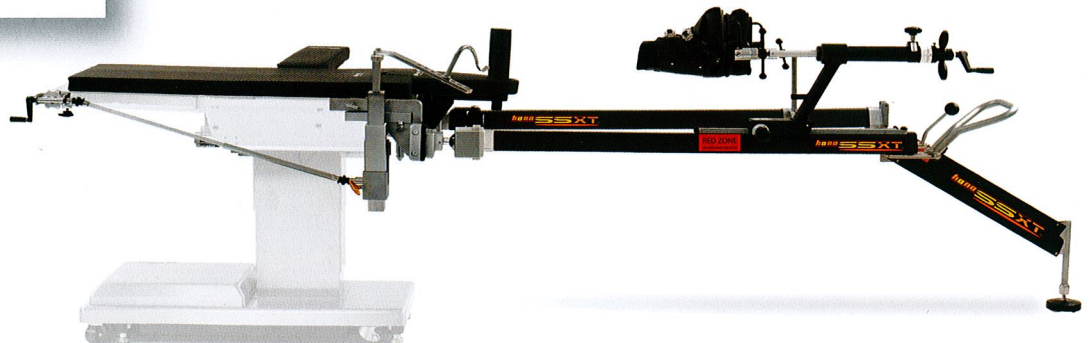

*Additional fit verification may be required for your specific general surgery table.

ORDERING INFORMATION

hanaSSXT™ For Total Hip Arthroplasty
and Hip Arthroscopy.....6870

6870 STANDARD COMPONENTS

Perineal Post, Large Diameter.....3369
Accessory Clamp..... 5855-61
Traction Boots, Small Left & Right..... 6850-485
Traction Boots, Large Left & Right..... 6850-486
Spar Mount Assembly, SSXT..... 6870-10
Upright Assembly, Well Leg..... 6870-29
Table Pad Body, SSXT..... 6870-39
Table Pad Perineal, SSXT..... 6870-40
Well Leg Spar Assembly..... 6870-6
Traction Spar Assembly..... 6870-7
Traction Hook Extension¹..... 6875-399
Perineal Post, Adult..... 6850-413
Traction Lifter..... 6850-480
Femoral Hook, Right¹..... 6850-918
Femoral Hook, Left¹..... 6850-919
Hand Crank Assembly..... 6870-41
Saddle Assembly Traction Lifter..... 6870-45
Femoral Hook Support¹, SSXT..... 6870-90

SPECIFICATIONS

Table Extension Length..... 87in. (221cm)
Average Overall Length.....130in. (330cm)
Operative Leg Spar Articulation
.....28° up-35° down
.....20° adduction-45° abduction
Patient Weight Capacity.....350lbs.(159.1kgs)

MIZUHOSI®

MIZUHO OSI
30031 AHERN AVENUE
UNION CITY, CA 94587/USA
(510) 429-1500
(800) 777-4674
FAX (888) 329-4674
FAX (510) 429-8500

MIZUHOSI.COM
NEWHIPNEWS.COM

ISO certified by TÜV Rheinland® of North America



Note: Mizuho OSI is constantly improving its products. All specifications are subject to change without notice. ¹ Patent US # 7,824,353 B2

Mizuho OSI is a Delaware Corporation. Manufactured in the USA

©Mizuho OSI LIT NW0590 Rev A 12/10

ORDERING INFORMATION

hanaSSXT™ Hip Arthroscopy Unit.....6874

6874 STANDARD COMPONENTS

Large Diameter Perineal Post..... 3369
Accessory Clamp..... 5855-61
Traction Boots, Small Left & Right..... 6850-485
Traction Boots, Large Left & Right..... 6850-486
Spar Mount Assembly, SSXT..... 6870-10
Upright Assembly Well Leg..... 6870-29
Table Pad Body, SSXT..... 6870-39
Table Pad Perineal, SSXT..... 6870-40
Well Leg Spar Assembly..... 6870-6
Traction Spar Assembly..... 6870-7
Traction Hook Extension..... 6875-399

OPTIONAL ACCESSORIES

hana® Cart..... 6875-500
Leg Support, Accessory..... 6850-170
Accessory Clamp..... 5855-61
Traction Boots X-Large Left & Right..... 6850-487
Femur Lift Accessory Kit (for 6874 only)..... 6870-1